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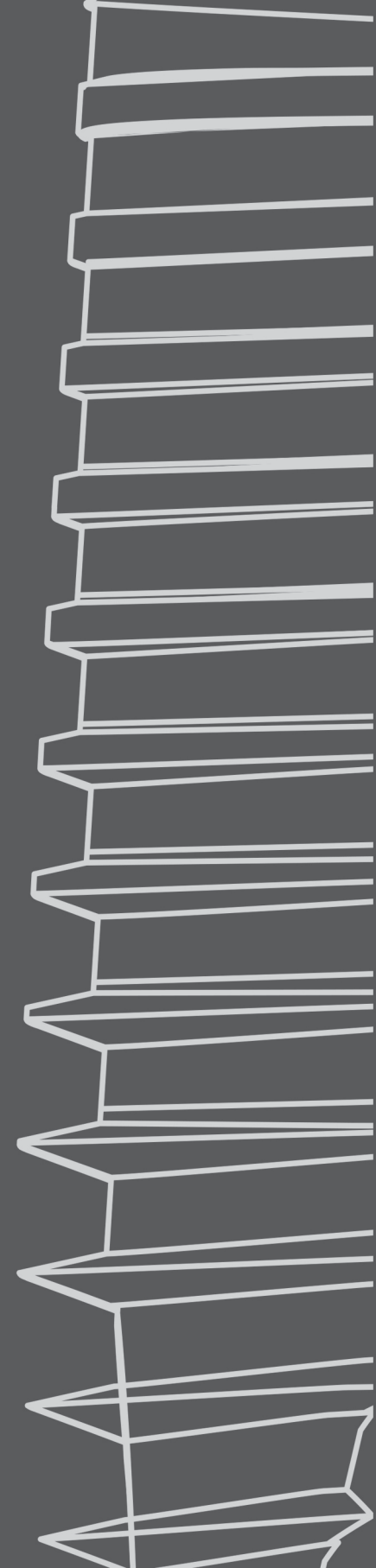
Studied by Buser (2004) , Zhao et Schwartz (2005) , Carl E Misch

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Studied by Huzeler M, Fickl S, Zuhr O , wachtel HC in J Oral Maxilofac Surg 2008 Oct 66

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
Diameter of Apex : Ø 2,1 mm


## Implant : ID<sup>®</sup>CAM

Reduces the risk of perforating the  
Schneiderian membrane.
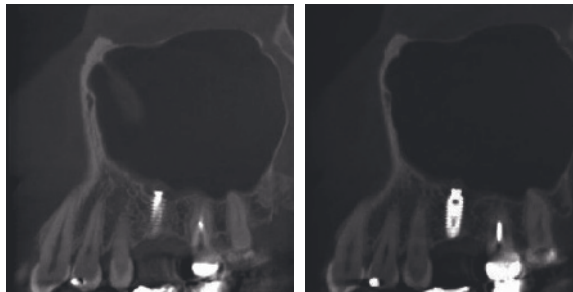
## Others Implants

The implant perforated the membrane.

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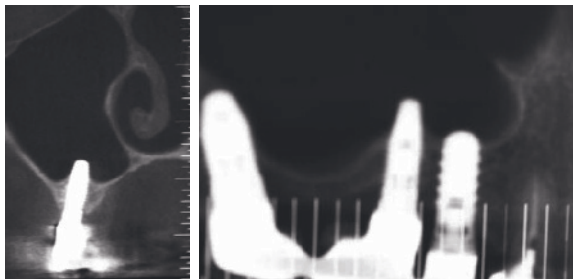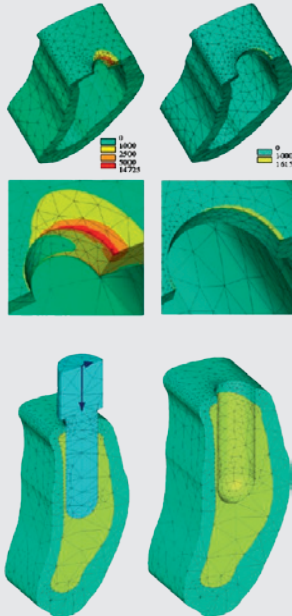
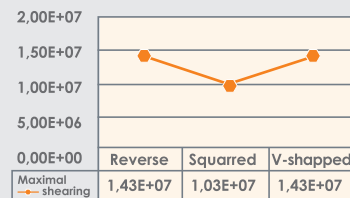
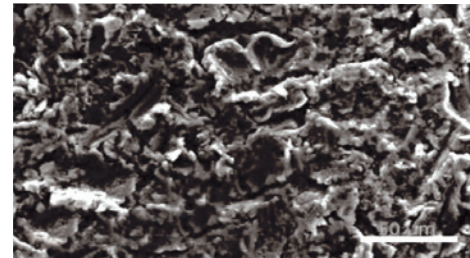
SMA®

Sandblasted-Macro Grits Acid Etched  
Neutralized by titanium Oxide

## Bone density for IDCAM ST


Type I

Type II

Bone quality classification of  
Lekholm & Zarb.

## Bone density for IDCAM M

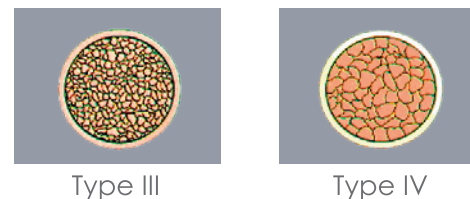Type III

Type IV

Bone quality classification of  
Lekholm & Zarb.

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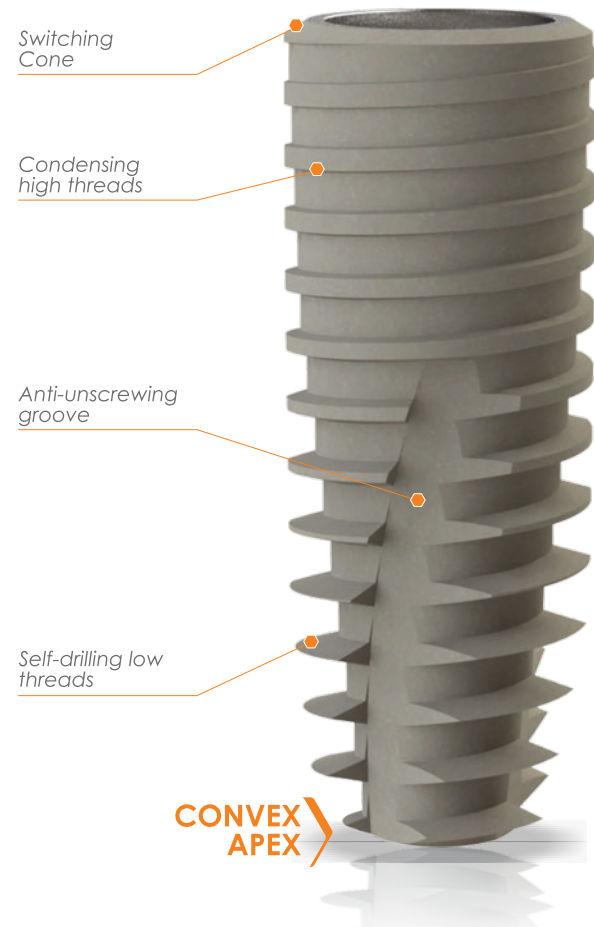
It is recommended to set as many implants as lacking natural roots in the patient mouth in order to secure the lasting of the prosthesis.

The length choice and diameter of implant must be based on the bone density determined by the CT- scan.

## IDCAM ST IMPLANT RANGE





| IMPLANT TYPE | LENGTH<br>Color code* | DIAMETER | PLATFORM | P/N       |
|--------------|-----------------------|----------|----------|-----------|
| STANDARD     | 8 mm ●                | 3,7 mm   | 3,6 mm   | IDCS0835  |
| STANDARD     | 8 mm ●                | 4,2 mm   | 3,6 mm   | IDCS0842  |
| STANDARD     | 8 mm ●                | 5,2 mm   | 3,6 mm   | IDCS0852  |
| STANDARD     | 10 mm ●               | 3,7 mm   | 3,6 mm   | IDCST1035 |
| STANDARD     | 10 mm ●               | 4,2 mm   | 3,6 mm   | IDCST1042 |
| STANDARD     | 10 mm ●               | 5,2 mm   | 3,6 mm   | IDCST1052 |
| STANDARD     | 12 mm ●               | 3,7 mm   | 3,6 mm   | IDCST1235 |
| STANDARD     | 12 mm ●               | 4,2 mm   | 3,6 mm   | IDCST1242 |
| STANDARD     | 12 mm ●               | 5,2 mm   | 3,6 mm   | IDCST1252 |
| STANDARD     | 15 mm ●               | 3,7 mm   | 3,6 mm   | IDCST1535 |
| STANDARD     | 15 mm ●               | 4,2 mm   | 3,6 mm   | IDCST1542 |
| STANDARD     | 15 mm ●               | 5,2 mm   | 3,6 mm   | IDCST1552 |

\*On each implant's packaging there is a small colored sticker to match with the implant height. The code for each color is related to the one found on the RBS conical and TURBOdrill® drills for the IDCAM implants.

- Length: 8mm
- Length: 10mm
- Length: 12mm
- Length: 15mm

Morse taper and identical connection for all IDCAM range.

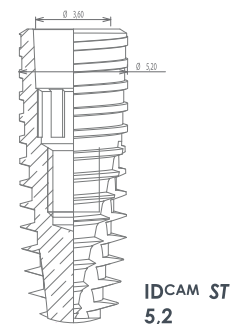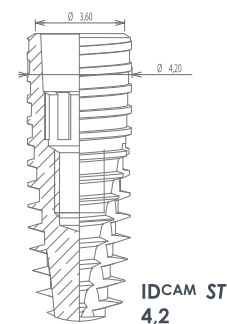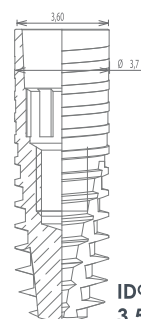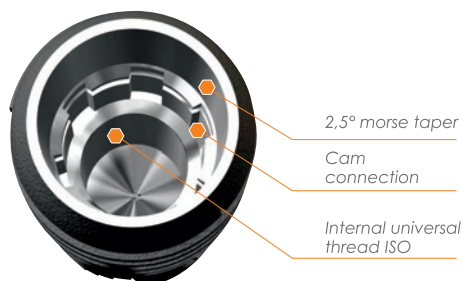SURGICAL  
SETS  
&  
DRILLS

**CONTENTS :** ®  
12 TURBOdrill  
4 RBS Ø2mm drill

1 pointing drill  
(L15 Ø1,8  
with depth marks  
P/N : 1518PT)

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

Hex-tipped screwdriver,  
short  
P/N 0046

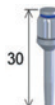Hex-tipped screwdriver,  
long  
P/N 0146

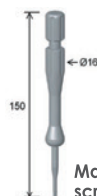Manual screwdriver to  
screw the implant  
P/N 0846

Screwdriver for relating  
instruments (short model)  
P/N 0014

Screwdriver for relating  
instruments (long  
model) P/N 0114

### Instruments

Instruments with  
dental shank, square-  
tipped instrument  
(short) P/N 1014

Instruments with  
dental shank, square-  
tipped instrument  
(long) P/N 1114

Instruments with  
dental shank, hex-  
tipped instrument  
(short) P/N 1046

Instruments with  
dental shank, hex-  
tipped instrument  
(long) P/N 1146

Manual thread tap,  
hardened stainless  
steel P/N TAR M2


Instrument  
extension  
P/N 406


**Double Use :**  
- Abutment remover for  
tapered implants (ID<sup>CAM</sup>  
and ID<sup>ALL</sup> product lines)  
- New square-tipped for  
for relating instruments  
(as P/N 0114)

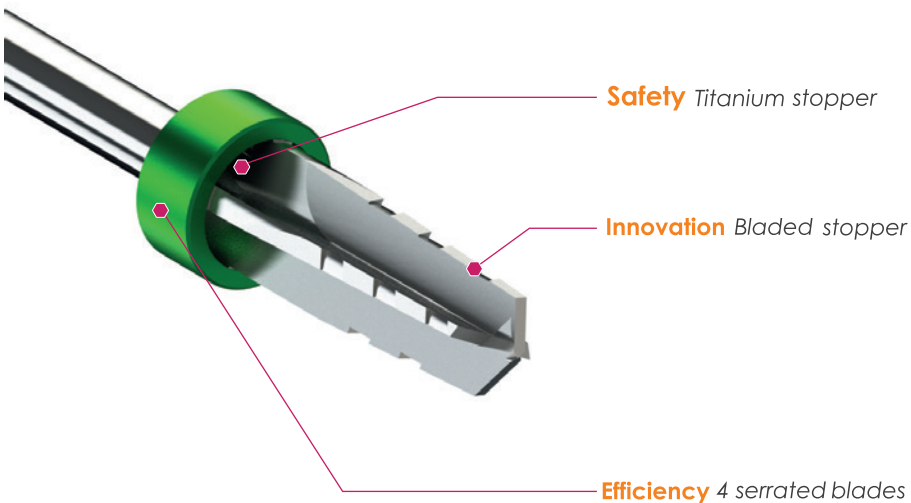ABUNDANT  
IRRIGATION DOWN  
INTO THE SOCKET

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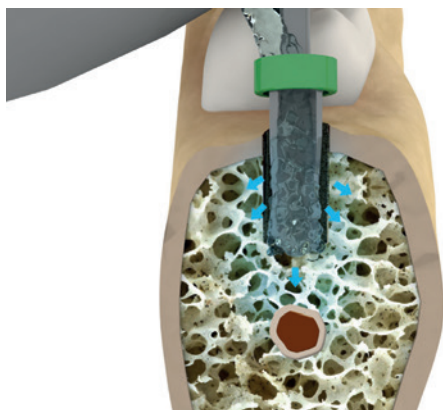
Press


Remove

➔ Pick up the closing cap from the packaging.

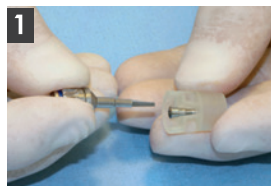
Take the packaging

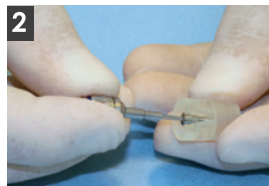
Pick

➔ Pick up the healing cap from the lower part of the packaging.

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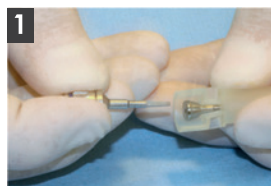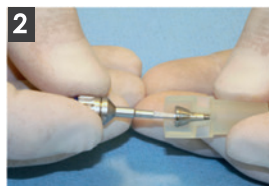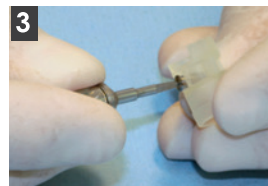
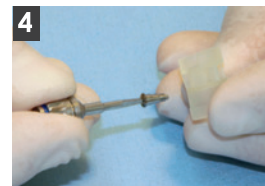
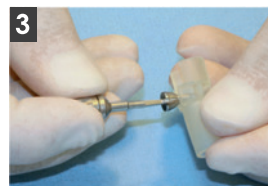
Remove

### USEFUL TIPS

1. Clean the instruments in an ultrasonic bath during 15 minutes for decontamination.
2. Rinse abundantly.
3. Dry perfectly.
4. Sterilize (dry heat) at 150° C during 90 minutes or at the autoclave at 135° during 20 minutes.

**Immediate Implant Placement Into Fresh Extraction Sites Using Single-Drilling Bur and Two Loading Procedures: Follow-Up Results**

Raphaël Bettach, DDS, Silvio Taschieri, MD, DDS, Carmen Mortellaro, MD, DDS, and Massimo Del Fabbro, MSc, PhD  
*The Journal of Craniofacial Surgery- Volume 29, Number 8, November 2018*

**Osteotomy at Low-Speed Drilling without Irrigation Versus High-Speed Drilling with Irrigation: an Experimental Study**

João GASPARG, Gonçalo BORRECHO Pedro OLIVEIRA, Francisco SALVADO, José MARTINS dos SANTOS  
*Acta Med Port 2013 May-Jun;26(3):231-236*

**Heat generation and drill wear during dental implant site preparation: systematic review**

S.C. Möhlhenrich, A. Modabbera, T. Steiner, D.A. Mitchell, F. Hölzle  
*Br J Oral Maxillofac Surg (2015)*

**Evaluation of the insertion torque, implant stability quotient and drilled hole quality for different drill design: an in vitro Investigation**

Sergio Alexandre Gehrke, José Luis Calvo Guirado, Raphaël Bettach, Massimo Del Fabbro, Carlos Perez-Albacete Martinez, Jamil Awad Shibli  
*Clin. Oral Impl. Res. 00, 2016, 1-7*

**Clinical evaluation of submerged and non-submerged implants for posterior single-tooth replacements: a randomized split-mouth clinical trial**

S. K. Nemli, M. B. Güngör, C. Aydın, H. Yılmaz, I. Türkcan, H. Demirköprülü - Department of Prosthodontics, Faculty of Dentistry, Gazi University, Ankara, Turkey  
*Int. J. Oral Maxillofac. Surg. 2014; 43: 1484-1492*

**Clinical and radiographic evaluation of new dental implant system: Results of a 3 year prospective study**

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**Implant Survival after Preparation of the Implant Site Using a Single Bur: A Case Series**

Raphaël Bettach, DDS; Silvio Taschieri, MD, DDS; Gilles Boukhris, DDS; Massimo Del Fabbro, BSc, PhD  
*Clinical Implant Dentistry and Related Research, (2013)*

**Temperature Changes in Cortical Bone after Implant Site Preparation Using a Single Bur versus Multiple Drilling Steps: An In Vitro Investigation**

Sergio Alexandre Gehrke, DDS; Raphaël Bettach, DDS; Silvio Taschieri, MD, DDS; Gilles Boukhris, DDS; Stefano Corbelli, DDS, PhD; Massimo Del Fabbro, BSc, PhD  
*Clinical Implant Dentistry and Related Research, (2013)*

**Clinical Performance of a Highly Porous Beta-TCP as the Grafting Material for Maxillary Sinus Augmentation**

Raphael Bettach, DDS, Bernard Guillaume, MD, Silvio Taschieri, MD, DDS, and Massimo Del Fabbro, BSc, PhD  
*Implant Dentistry, (2014)*

Find all the publications on :  
[www.idi-dental.com/documentation](http://www.idi-dental.com/documentation)

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Open Tray (Anti Rotational): 2004  
Close Tray (Anti Rotational): 2004F  
Closed Tray (Rotational): 0221

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

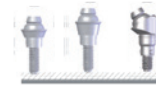
3- The **FUNCTIONAL ASPECT IS PRESERVED** thanks to the **asymmetric shoulder shape**, including a lower vestibular zone than the proximal zone (papilla).
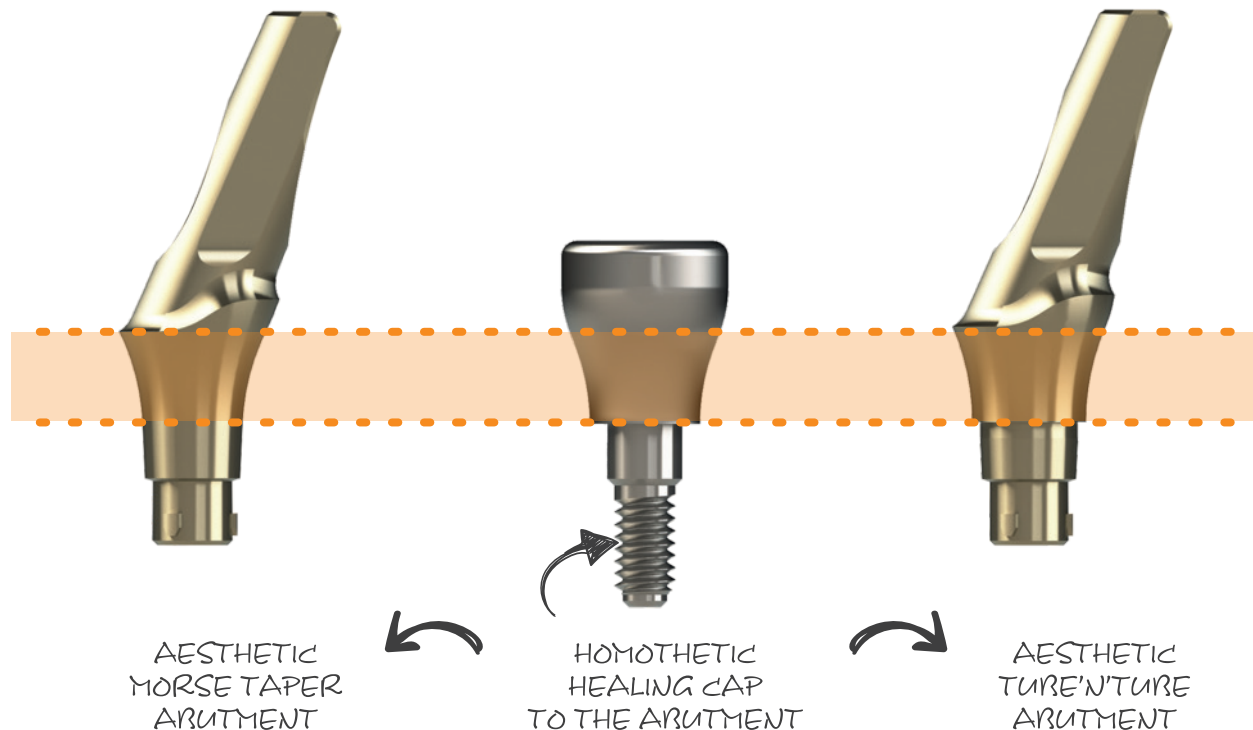
| IMPLANT TYPE | LENGTH<br>Color code* | DIAMETER      | PLATFORM | P/N       |
|--------------|-----------------------|---------------|----------|-----------|
| STANDARD     | 8 mm ●                | <b>3,7 mm</b> | 3,6 mm   | IDCS0835  |
| STANDARD     | 8 mm ●                | <b>4,2 mm</b> | 3,6 mm   | IDCS0842  |
| STANDARD     | 8 mm ●                | <b>5,2 mm</b> | 3,6 mm   | IDCS0852  |
| STANDARD     | 10 mm ●               | <b>3,7 mm</b> | 3,6 mm   | IDCST1035 |
| STANDARD     | 10 mm ●               | <b>4,2 mm</b> | 3,6 mm   | IDCST1042 |
| STANDARD     | 10 mm ●               | <b>5,2 mm</b> | 3,6 mm   | IDCST1052 |
| STANDARD     | 12 mm ●               | <b>3,7 mm</b> | 3,6 mm   | IDCST1235 |
| STANDARD     | 12 mm ●               | <b>4,2 mm</b> | 3,6 mm   | IDCST1242 |
| STANDARD     | 12 mm ●               | <b>5,2 mm</b> | 3,6 mm   | IDCST1252 |
| STANDARD     | 15 mm ●               | <b>3,7 mm</b> | 3,6 mm   | IDCST1535 |
| STANDARD     | 15 mm ●               | <b>4,2 mm</b> | 3,6 mm   | IDCST1542 |
| STANDARD     | 15 mm ●               | <b>5,2 mm</b> | 3,6 mm   | IDCST1552 |

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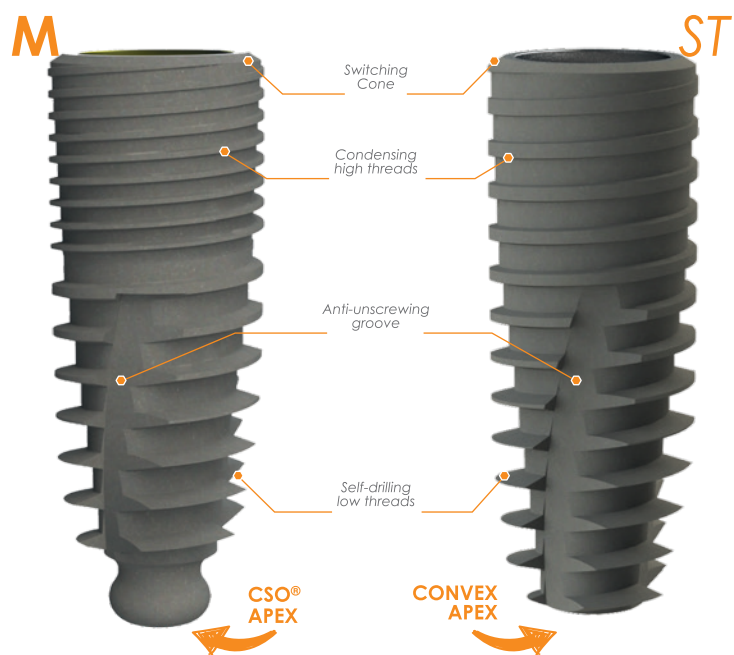
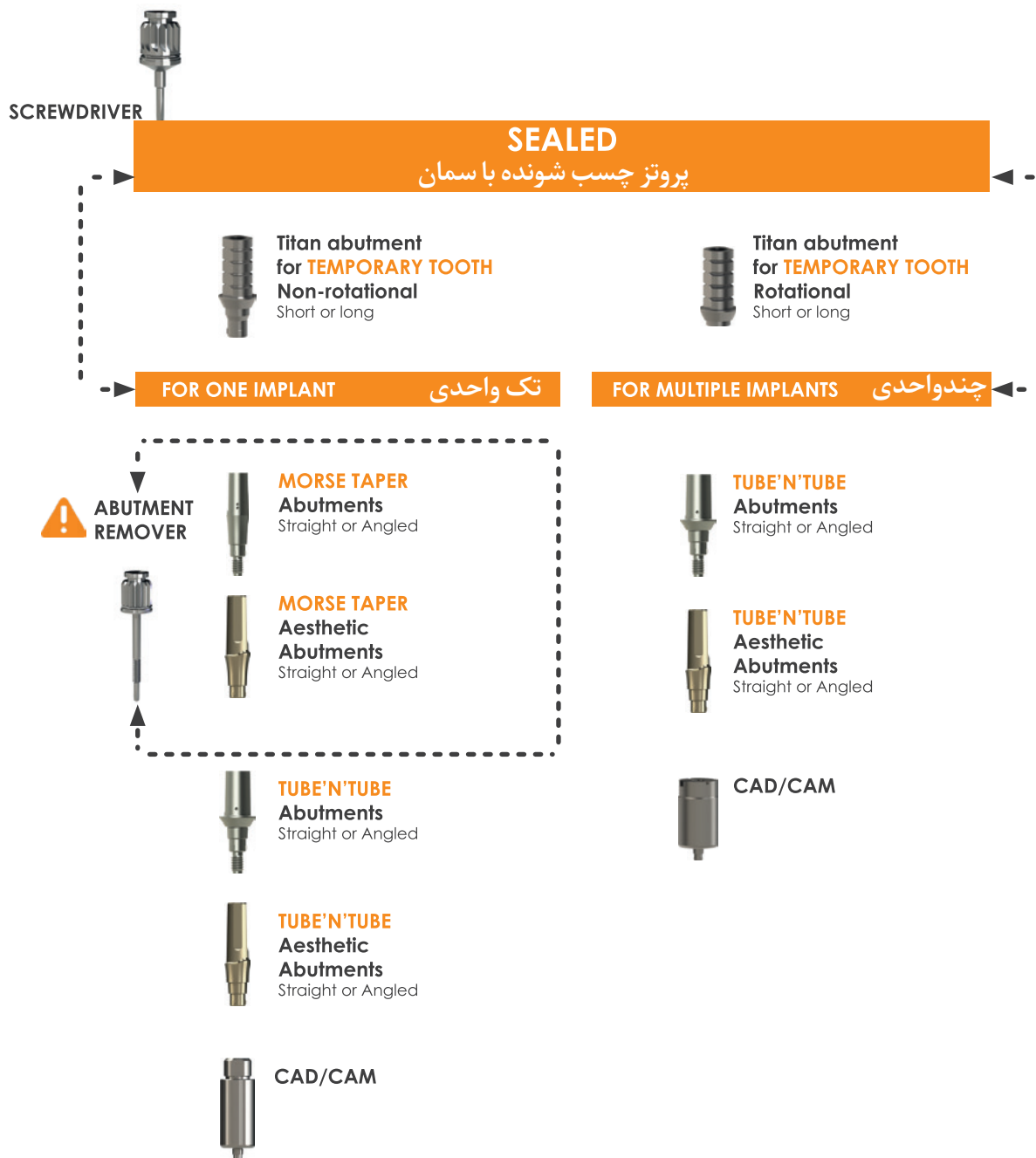
2 PROSTHETIC CONNEXION TYPE.  
AN ANSWER TO EACH CASE.

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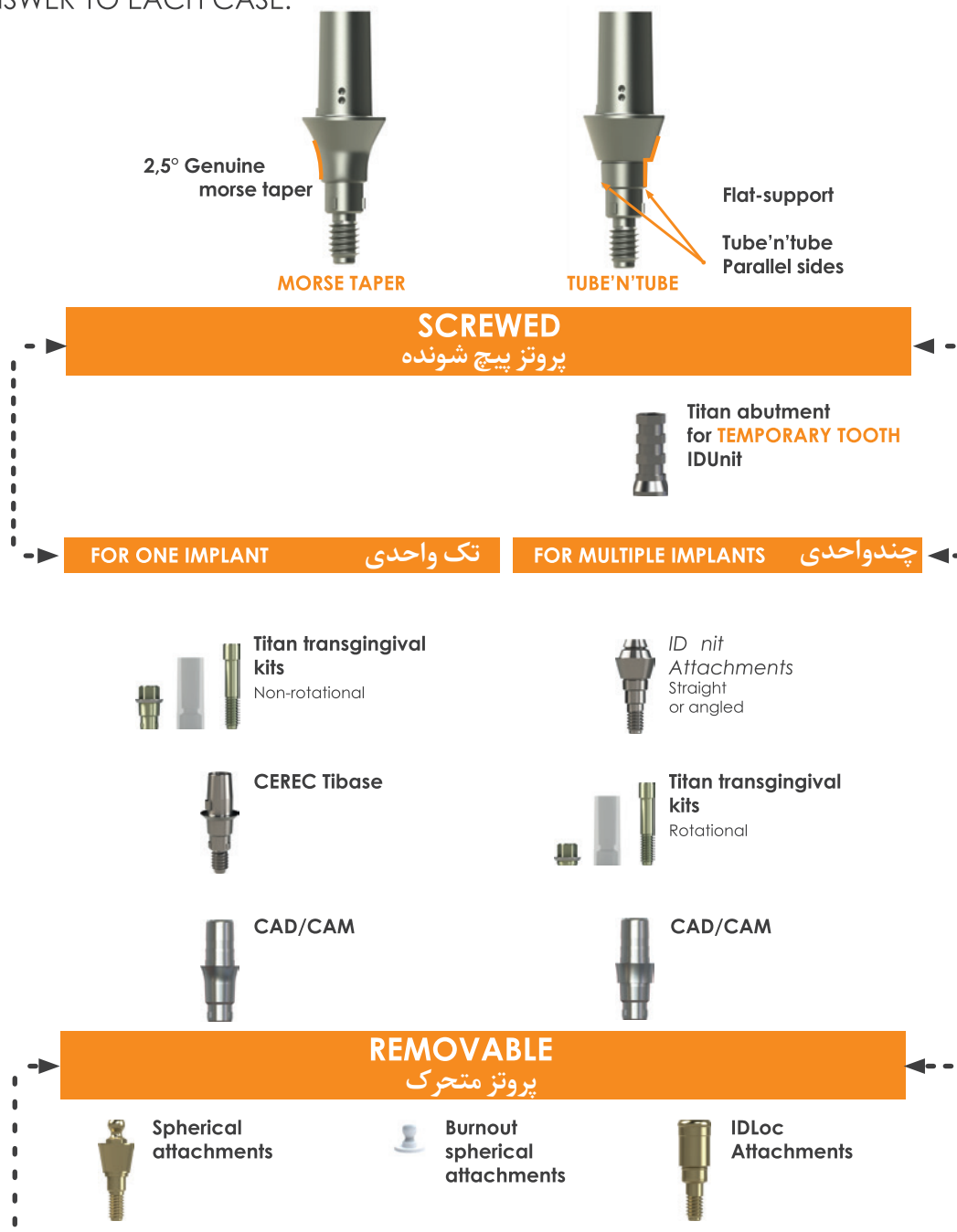## → PROSTHETIC COMPONENTS

| IMPRESSION CYLINDER <b>NON-ROTATIONAL</b>                                         |                                    |        |
|-----------------------------------------------------------------------------------|------------------------------------|--------|
|  | Pick-up technique                  | 2004   |
|  | Closed tray technique              | 2004F  |
|  | Closed tray technique, long        | 2004FL |
|  | Pick-up technique, narrow          | 2004N  |
|  | Closed tray technique, narrow long | 2004NL |
|  | Pick-up technique <b>Plastic</b>   | 2004P  |

| IMPRESSION CYLINDER <b>ROTATIONAL</b>                                              |                                            |       |
|------------------------------------------------------------------------------------|--------------------------------------------|-------|
|  | Conical                                    | 0220C |
|  | Rotational                                 | 0221  |
|  | Rotational long, special post-extractional | 0221L |

| RETAINING SCREWS                                                                   |                                                                                   |      |
|------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|------|
|  | Retaining screw for tapped screw-retained elements<br><b>25N.cm maximum</b>       | 0211 |
|  | Retaining screw for ID <sup>CAM</sup> and ID <sup>ALL</sup> <b>25N.cm maximum</b> | 0214 |
|  | Long head screw for laboratory for all ranges, <b>25N.cm maximum</b>              | 0217 |
|  | Short head, <b>25N.cm maximum</b>                                                 | 0219 |

| SCANBODY                                                                            |                                              |         |
|-------------------------------------------------------------------------------------|----------------------------------------------|---------|
|   | Scanbody CEREC® Bluecam <b>S type</b> (by 5) | 6431295 |
|  | Scanbody CEREC® Bluecam <b>L type</b> (by 5) | 6431303 |
|  | Scanbody CEREC® Omnicam <b>S type</b> (by 5) | 6431311 |
|  | Scanbody CEREC® Omnicam <b>L type</b> (by 5) | 6431329 |

| IMPLANT ANALOG                                                                      |                           |      |
|-------------------------------------------------------------------------------------|---------------------------|------|
|  | Implant analog<br>Ø3,6 mm | 0223 |

| CEREC TIBASE                                                                         |                             |      |
|--------------------------------------------------------------------------------------|-----------------------------|------|
|  | Titan base for CEREC Ø3,6mm | 7336 |
|  | Titan base for CEREC Ø3,1mm | 7436 |

## → THE TEMPORARY PROSTHESIS

| TITAN ABUTMENTS FOR A TEMPORARY TOOTH<br>FOR ONE IMPLANT RESTORATION                |                                                                           |       |
|-------------------------------------------------------------------------------------|---------------------------------------------------------------------------|-------|
|  | Non-rotational - Ø3,6 mm                                                  | 0206  |
|  | Temporary titan abutment, long (post-extractional for switching platform) | 0206L |

| TITAN ABUTMENTS FOR A TEMPORARY TOOTH<br>FOR MULTIPLE IMPLANTS RESTORATION           |                                                                           |       |
|--------------------------------------------------------------------------------------|---------------------------------------------------------------------------|-------|
|  | Rotational - Ø3,6 mm                                                      | 0208  |
|  | Temporary titan abutment, long (post-extractional for switching platform) | 0208L |

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\*H.T. : Transgingival height

**USE THE MORSE TAPER ABUTMENT REMOVER,**  
if needed.  
p107

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USE THE MORSE TAPER  
ABUTMENT REMOVER,  
if needed.

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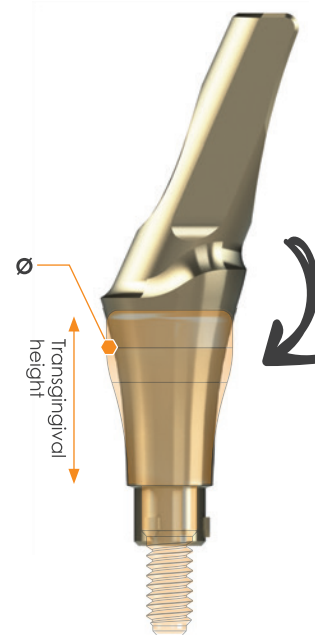
### BURNOUT CYLINDER

|                                                                                   |                                                                    |        |
|-----------------------------------------------------------------------------------|--------------------------------------------------------------------|--------|
|  | Non-rotational with shoulder in nylon Ø4,8 mm, for temporary tooth | 022602 |
|-----------------------------------------------------------------------------------|--------------------------------------------------------------------|--------|

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#### IDUNIT IMPRESSION COPYINGS

|                                                                                   |                         |      |
|-----------------------------------------------------------------------------------|-------------------------|------|
|  | Closed tray technique   | 321  |
|  | Pick-up technique       | 322  |
|  | Pick-up technique, long | 322L |

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| THE BURNOUT ELEMENTS                                                              |                                                     |      |
|-----------------------------------------------------------------------------------|-----------------------------------------------------|------|
|  | Spherical attachments                               | 9222 |
|  | Paralleling guide for burnout spherical attachments | 9223 |
|  | Burnout connector bar <b>by 3</b>                   | 0931 |
|  | Nylon cap                                           | 0025 |


### IDLOC SYSTEM

| IDLOC IMPLANT ANALOG                                                                |                |     |
|-------------------------------------------------------------------------------------|----------------|-----|
|  | Implant analog | 433 |

| IDLOC IMPRESSION COPYING                                                            |         |     |
|-------------------------------------------------------------------------------------|---------|-----|
|  | Plastic | 432 |

| IDLOC ATTACHMENT BOX                                                                |                                                 |        |
|-------------------------------------------------------------------------------------|-------------------------------------------------|--------|
|  | Female part (Locator)<br><b>Thickness 2,5mm</b> | LOCFEM |

| BOXES FOR SPHERICAL ATTACHMENTS                                                    |                                                                   |        |
|------------------------------------------------------------------------------------|-------------------------------------------------------------------|--------|
|  | O'ring <b>height:3,5mm</b><br><b>Ø external 5mm</b>               | 0122   |
|  | Nitrile retaining ring soft (white)                               | 0120NB |
|  | O'ring retaining ring for O'ring attachment <i>medium</i> (red)   | 0120SR |
|  | O'ring retaining ring for O'ring attachment <i>strong</i> (black) | 0120NN |
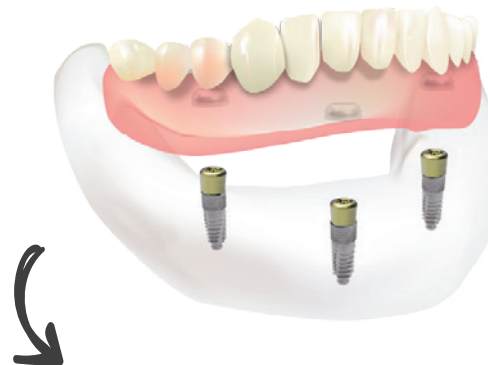
| IDLOC ATTACHMENTS                                                                     |                                    |       |
|---------------------------------------------------------------------------------------|------------------------------------|-------|
|  | Transgingival height <b>1 mm</b>   | L3601 |
|  | Transgingival height <b>2,5 mm</b> | L3602 |
|  | Transgingival height <b>4 mm</b>   | L3604 |
|  | Transgingival height <b>6 mm</b>   | L3606 |

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### TITANIUM SCANBODY


|                                                                                   |                                              |      |
|-----------------------------------------------------------------------------------|----------------------------------------------|------|
|  | CAD/CAM impression copying<br>Cam connection | SBCT |
|-----------------------------------------------------------------------------------|----------------------------------------------|------|

### PREMILLED ABUTMENT Ø12

|                                                                                   |             |       |
|-----------------------------------------------------------------------------------|-------------|-------|
|  | Conical     | PMC12 |
|  | Tube'n'Tube | PMP12 |

### PREMILLED ABUTMENT Medentika Type Ø12

|                                                                                    |             |        |
|------------------------------------------------------------------------------------|-------------|--------|
|   | Conical     | PMC12M |
|  | Tube'n'Tube | PMP12M |
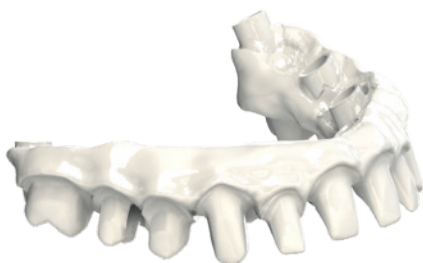
### IDUNIT

|                                                                                     |                                      |      |
|-------------------------------------------------------------------------------------|--------------------------------------|------|
|  | IDUnit Tibase + screw P/N 0216       | TU   |
|  | IDUnit CAD/CAM impression<br>copying | 321P |

### TIBASES

#### TUBE'N'TUBE ROTATIONAL CONNECTION Ø 3.6MM

|                                                                                   |                                      |       |
|-----------------------------------------------------------------------------------|--------------------------------------|-------|
|  | Height 0,5 mm<br>With screw P/N 0211 | TR360 |
|  | Height 2 mm<br>With screw P/N 0211   | TR362 |
|  | Height 4 mm<br>With screw P/N 0211   | TR364 |

#### CONICAL CONNECTION Ø 3.6MM

|                                                                                   |                                      |       |
|-----------------------------------------------------------------------------------|--------------------------------------|-------|
|  | Height 0,5 mm<br>With screw P/N 0211 | TC320 |
|  | Height 2 mm<br>With screw P/N 0211   | TC322 |
|  | Height 4 mm<br>With screw P/N 0211   | TC324 |

#### TUBE'N'TUBE CAM CONNECTION Ø 3.6MM

|                                                                                     |                                      |       |
|-------------------------------------------------------------------------------------|--------------------------------------|-------|
|  | Height 0,5 mm<br>With screw P/N 0211 | TP360 |
|  | Height 2 mm<br>With screw P/N 0211   | TP362 |
|  | Height 4 mm<br>With screw P/N 0211   | TP364 |

IMPLANTS  
DIFFUSION INTERNATIONAL  
Since 1987

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
Platform  
Switch

Double twist  
threads

Penetrating and  
condensing  
threads

Anti-unscrewing  
groove

**CATCH  
BASE**

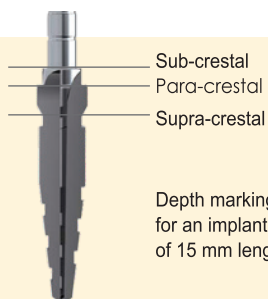Depth markings  
for an implant  
of 15 mm length.

There is one drill for each implant length.

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
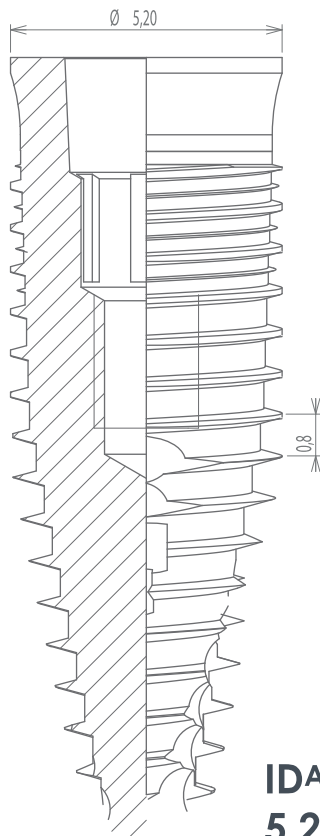
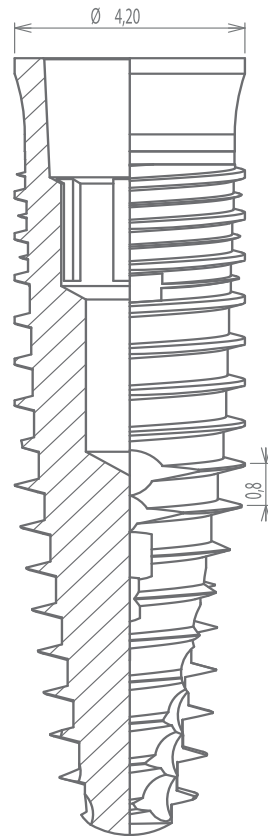
All the dimensions are in millimeters.

# SURGICAL DRILL SET & DRILLS FOR THE ID<sup>ALL</sup> IMPLANTS

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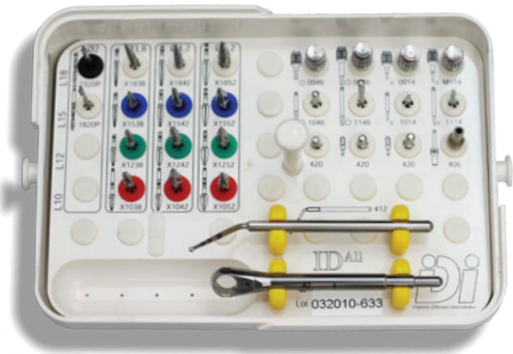
- CONTENTS:**
- 2 pilot drills
  - 12 cylindro-tapered drills with 4 bladed edges
  - 11 relating instruments

**Markings**  
Markings on the shank indicate the diameter of the drill.

### RANGE OF ID<sup>ALL</sup> DRILLS

|          | DRILLING LENGTH<br>Color code* | P/N   |
|----------|--------------------------------|-------|
| Ø 3,8 mm | 10                             | X1038 |
| Ø 4,2 mm | 10                             | X1042 |
| Ø 5,2 mm | 10                             | X1052 |
| Ø 3,8 mm | 12                             | X1238 |
| Ø 4,2 mm | 12                             | X1242 |
| Ø 5,2 mm | 12                             | X1252 |
| Ø 3,8 mm | 15                             | X1538 |
| Ø 4,2 mm | 15                             | X1542 |
| Ø 5,2 mm | 15                             | X1552 |
| Ø 3,8 mm | 18                             | X1838 |
| Ø 4,2 mm | 18                             | X1842 |
| Ø 5,2 mm | 18                             | X1852 |
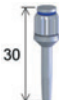
|                                                                                                                                                                                 |                                                                                                                                                                                    |                                                                                                                                                                                 |                                                                                                                                                                                   |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  <p><b>Ratchet</b></p> <p>P/N 415</p>                                                         |  <p><b>Instrument extension</b></p> <p>P/N 406</p>                                               |  <p><b>Depth tester</b></p> <p>P/N 412</p>                                                    |  <p><b>Implant paralleling help by 3</b></p> <p>P/N 420</p>                                   |
|  <p><b>Handpiece screwdriver for prosthetics, length: 30 mm</b></p> <p>P/N 0114</p>          |  <p><b>Contra angle screwdriver for prosthetics, length: 30 mm</b></p> <p>P/N 1114</p>          |  <p><b>Handpiece screwdriver for screwing the implant, length: 22 mm</b></p> <p>P/N 0014</p> |  <p><b>Contra angle screwdriver for prosthetics, length: 22 mm</b></p> <p>P/N 1014</p>       |
|  <p><b>Handpiece screwdriver for screwing the implant, length: 30 mm</b></p> <p>P/N 0146</p> |  <p><b>Contra angle screwdriver for screwing the implant, length: 30 mm</b></p> <p>P/N 1146</p> |  <p><b>Handpiece screwdriver for screwing the implant, length: 22 mm</b></p> <p>P/N 0046</p> |  <p><b>Handpiece screwdriver for screwing the implant, length: 32 mm</b></p> <p>P/N 1046</p> |

The ID<sup>ALL</sup> drills were developed and tested in several hospitals. They ease the implant placement of the ID<sup>ALL</sup> and ID<sup>ALL</sup>S ranges.

\*On each implant's packaging there is a small colored sticker to match with the implant height. The code for each color is related to the one found on the ID<sup>ALL</sup> drills for the ID<sup>ALL</sup> implants.

- Length: 10 mm
- Length: 12 mm
- Length: 15 mm
- Length: 18 mm

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### REMARKS

For a total edentation and in order to stabilize a removable denture, it is recommended to set a minimum of 4 ID<sup>SLIM</sup> implants at the mandible and 6 ID<sup>SLIM</sup> implants at the maxillary, as parallel as possible.

# THE RANGE ID<sup>SLIM</sup>

## ID<sup>SLIM</sup> surgical set containing :

|  |                                                                               |       |
|--|-------------------------------------------------------------------------------|-------|
|  | Ratchet                                                                       | 415   |
|  | Instrument extension with dental shank                                        | 406   |
|  | Paralleling guide Ø2mm / Ø2,35 mm x 5                                         | 410   |
|  | Manual screwdriver to screw the ID <sup>SLIM</sup> implant Ø8mm, Length: 22mm | 0246  |
|  | Manual screwdriver to screw the ID <sup>SLIM</sup> implant Ø8mm, Length: 30mm | 0346  |
|  | Instrument with dental shank to screw the ID <sup>SLIM</sup> implant, short   | 1246  |
|  | Cylindrical drill L18 Ø2mm with depth marking                                 | 1820P |
|  | Cylindrical drill L18 Ø2,5mm with depth marking                               | 1825P |
|  | Cylindrical drill L18 Ø2,7mm with depth marking                               | 1827P |

## Implants


|  |                                     |           |
|--|-------------------------------------|-----------|
|  | Length 10mm - Height 1,5mm - Ø2,8mm | IDSLI1031 |
|  | Length 10mm - Height 3mm - Ø2,8mm   | IDSLI1033 |
|  | Length 12mm - Height 1,5mm - Ø2,8mm | IDSLI1231 |
|  | Length 12mm - Height 3mm - Ø2,8mm   | IDSLI1233 |
|  | Length 15mm - Height 1,5mm - Ø2,8mm | IDSLI1531 |
|  | Length 15mm - Height 3mm - Ø2,8mm   | IDSLI1533 |

## Prosthetic elements

|  |                                                             |      |
|--|-------------------------------------------------------------|------|
|  | Hexagonal abutment, straight                                | 3000 |
|  | Hexagonal abutment, 7° angled                               | 3007 |
|  | Hexagonal abutment, 15° angled                              | 3015 |
|  | Box for O'ring ball attachment Height: 3,5mm, external Ø5mm | 0122 |
|  | Retaining ring for O'ring attachment x3                     | 0120 |
|  | Impression copying                                          | 821  |
|  | Implant analog                                              | 823  |

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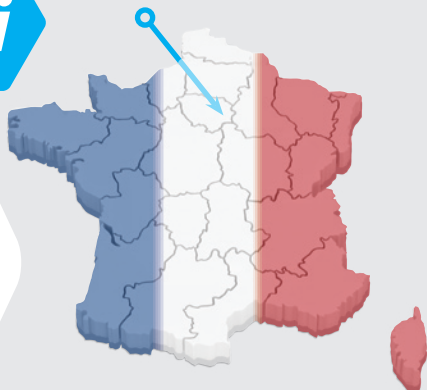idi offers a full Line of dental Implants meeting all your surgical requirements .

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
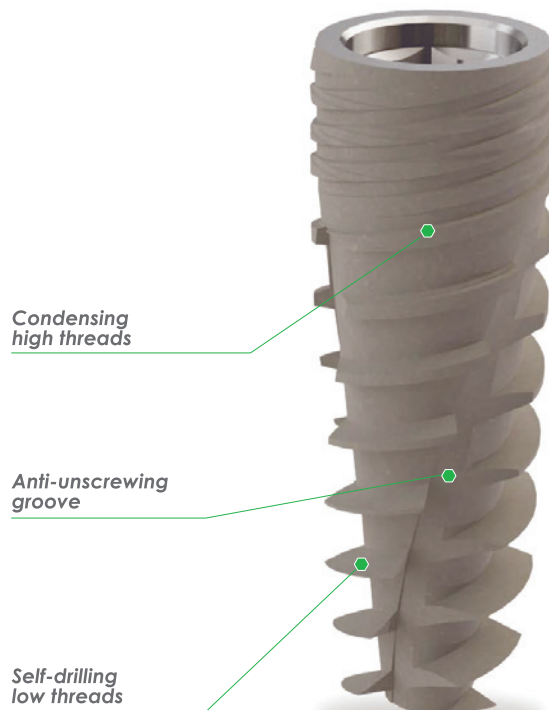
Anchorage  
of prosthetic  
stabilization

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

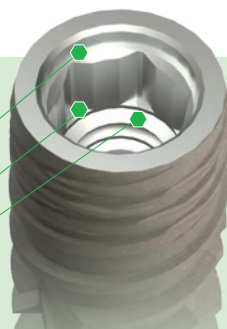
- Anti-unscrewing groove

- Progressive and condensing threads

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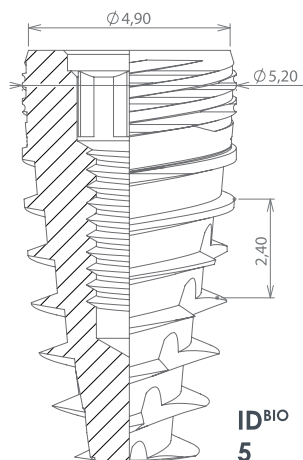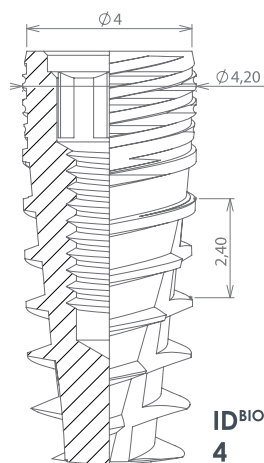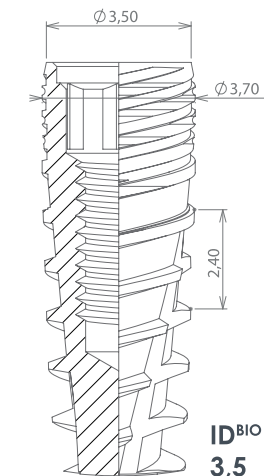### Golden tapped Hexagonal abutments

|  |                                                                       |       |       |       |
|--|-----------------------------------------------------------------------|-------|-------|-------|
|  | Straight Slim abutment for Platform switch + retaining screw P/N 0215 | 3400G | 3800G | 3800G |
|  | Abutment, straight + retaining screw P/N 0215                         | 3600G | 4000G | 5000G |
|  | Abutment, straight, with a 2 mm shoulder + retaining screw P/N 0215   | 3602G | 4002G | 5002G |
|  | Abutment, straight, with a 4 mm shoulder + retaining screw P/N 0215   | 3604G | 4004G | 5004G |
|  | Abutment, 15° angled + retaining screw P/N 0215                       | 3615G | 4015G | 5015G |
|  | Abutment, 23° angled + retaining screw P/N 0215                       | 3623G | 4023G | 5023G |
|  | Abutment, 0° to 18° + retaining screw P/N 0215                        | HTU1G | HTU2G | HTU5G |

### Retaining screws

|  |                                                                                  |      |
|--|----------------------------------------------------------------------------------|------|
|  | Golden Retaining screw for tapped screw-retained elements (Screw head: Ø 2,2 mm) | 0215 |
|  | For prosthesis: ≤ 25 N.cm maximum (Screw head: Ø 2,5 mm)                         | 1414 |
|  | For prosthesis: ≤ 25 N.cm maximum (short head Ø 2,5 mm)                          | 1413 |

### Burnout cylinders

|  |                                            |      |      |      |
|--|--------------------------------------------|------|------|------|
|  | Non rotational, + retaining screw P/N 1414 | 356S | 456S | 556S |
|  | Rotational, + retaining screw P/N 1414     | 358S | 458S | 558S |

### Titanium transgingival kits\*

|  |                                                                                     |       |       |       |
|--|-------------------------------------------------------------------------------------|-------|-------|-------|
|  | Non rotational element, 0,5 mm high Titanium basis & Burnout element Screw P/N 6140 | 6350H | 6400H | 6500H |
|  | Rotational element, 0,5 mm high Titanium basis & Burnout element Screw P/N 6140     | 6350R | 6400R | 6500R |
|  | Non rotational element, 1,5 mm high Titanium basis & Burnout element Screw P/N 6141 | 6351H | 6401H | 6501H |
|  | Rotational element, 1,5 mm high Titanium basis & Burnout element Screw P/N 6141     | 6351R | 6401R | 6501R |
|  | Non rotational element, 2,5 mm high Titanium basis & Burnout element Screw P/N 6142 | 6352H | 6402H | 6502H |
|  | Rotational element, 2,5 mm high Titanium basis & Burnout element Screw P/N 6142     | 6352R | 6402R | 6502R |

### CEREC

|  |                                                   |         |         |         |
|--|---------------------------------------------------|---------|---------|---------|
|  | Scanpost (small) + screw P/N 0211                 | 73HS    | 73HS    | 73HS    |
|  | Titanium abutment basis, for lab + screw P/N 0215 | 7340    | 7340    | 7350    |
|  | Omniscam Scanbody (small)                         | 6431311 | 6431311 | 6431311 |
|  | Bluecam Scanbody (small)                          | 6431295 | 6431295 | 6431295 |

### Zirconia TiBase

|  |                                          |      |      |      |
|--|------------------------------------------|------|------|------|
|  | Titanium abutment basis + screw P/N 0215 | 7535 | 7540 | 7550 |
|--|------------------------------------------|------|------|------|

### Gold cylinders

|  |                                                            |       |       |       |
|--|------------------------------------------------------------|-------|-------|-------|
|  | Cylinder, gold, non rotational to cast on + screw P/N 1414 | ASH35 | ASH40 | ASH50 |
|  | Cylinder, gold, rotational to cast on + screw P/N 1414     | ASR35 | ASR40 | ASR50 |

### IMPORTANT NOTICE:

ALL prosthetic elements have to be used with the screwdrivers P/N 0014, 0114, 1014, 1114 and 0148.

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| Boxes for spherical attachments                                                     |                                                |      |
|-------------------------------------------------------------------------------------|------------------------------------------------|------|
|    | O'ring, Height: 3,5 mm<br>External Ø 5 mm      | 0122 |
|   | O'ring retaining ring<br>for O'ring attachment | 0120 |
|  | Nylon box<br>for spherical attachment          | 0924 |

| Burnout spherical attachments                                                       |                                                       |      |
|-------------------------------------------------------------------------------------|-------------------------------------------------------|------|
|  | Burnout spherical attachment                          | 9222 |
|  | Paralleling guide<br>for burnout spherical attachment | 9223 |

| Connector bar                                                                       |                                 |      |
|-------------------------------------------------------------------------------------|---------------------------------|------|
|  | Burnout connector bar<br>(by 3) | 0931 |
|  | Nylon clip                      | 0025 |

| Magnetic screws                                                                     |                    |      |         |
|-------------------------------------------------------------------------------------|--------------------|------|---------|
|  | Screw, REDEIM type | 935  | 941 951 |
|  | Magnetic           | 0940 |         |

ITPGOLD

IMPLANTS DIFFUSION INTERNATIONAL  
Since 1987  
Made in FRANCE

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IDCAM ST

The company was founded in 1987 under the impulse of a team of maxillo-facial surgeons and dentists at the early stages of modern implantology. these pioneers developed the IDI implantology system in the Tenon hospital in Paris (France).

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
|                                                                                                                                                                                                                                                                                         |                                                                                                                                                                                                                                                                          |                                                                                                                                                                                                                                                    |                                                                                                                                                                                                                                                 |                                                                                                                                                                                                                                                 |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  <b>O.N.B</b><br> 2.7mm 15-18mm<br>3.2mm 8-10-12-15mm<br>3.7mm 8-10-12-15mm<br>5 mm 8-10-12-15mm<br>6 mm 8-10-12-15mm |  <b>ITPGOLD</b><br> 3.2mm 8-10-12-15mm<br>3.7mm 8-10-12-15mm<br>5 mm 8-10-12-15mm<br>6 mm 8-10-12-15mm |  <b>IDMAX</b><br> 4mm 6-8-10-12-15mm<br>5mm 6-8-10-12-15mm<br>6mm 6-8-10-12-15mm |  <b>IDOH</b><br> 4mm 8-10-12-15mm<br>5mm 8-10-12-15mm<br>6mm 8-10-12-15mm |  <b>IDOT</b><br> 4mm 8-10-12-15mm<br>5mm 8-10-12-15mm<br>6mm 8-10-12-15mm |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|

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|                                                                                                                                                                   |                                                                                                                                                               |                                                                                                                                                                                 |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  Closing cap<br>Order no. 341-441-541-641                                        |  Type REDEIM magnetic screw, 4mm<br>Order no. 941-951-961                    |  Titanium sealed abutment<br>Order no. HSU                                                   |
|  Healing cap, 5mm high<br>Order no. 313-413-513-613                              |  Retaining screw for screw on prosthetic units, long<br>Order no. 914        |  Titanium screw abutment<br>Order no. 921                                                    |
|  Impression coping, straight (for solid units)<br>Order no. 721                  |  Retaining screw for screw on prosthetic units, long<br>Order no. 914C       |  Temporary cylinder rotationally secure (for a temporary tooth)<br>Order no. APPH40-44-50-60 |
|  Impression coping, conical (for solid units)<br>Order no. 721CO                 |  Burnout connector bar, 5bar<br>Order no. 93.1                               |  Temporary cylinder rotational (for a temporary tooth)<br>Order no. APPR40-44-50-60          |
|  Impression coping, rotationally secure (for solid units)<br>Order no. H904    |  Ball anchor matrix<br>Order no. 924M                                      |  Abutment, 5 angled<br>Order no. AM40-44-50-60(05)                                         |
|  Impression coping, pick-up technique (rotationally secure)<br>Order no. H904L |  Implant analogue<br>Order no. 323-423-523-623                             |  Abutment, 15 angled (aesthetic)<br>Order no. AM40-44-50-60(15)                            |
|  Transgingival extension<br>Order no. 325-425-525-625                          |  Burnout cylinder, rotational (3 pieces)<br>Order no. 318-418-518-618      |  Abutment, 23 angled (aesthetic)<br>Order no. AM40-44-50-60(23)                            |
|  Titanium ball anchor, long<br>Order no. 922L                                  |  Burnout cylinder, rotationally (1 piece)<br>Order no. H326-426-526-626    |  Cylinder, Gold rotational (to cast on)<br>Order no. ASR40-44-50-60                        |
|  Ball anchor, short<br>Order no. 922C                                          |  Screw-on abutment (0-18)<br>Order no. HTU3-4-5-6<br>Order no. HTU3-4-5-6G |  Cylinder, Gold rotationally secure (to cast on)                                           |
|  Type REDEIM magnet, 4mm<br>Order no. 940                                      |                                                                                                                                                               |                                                                                                                                                                                 |

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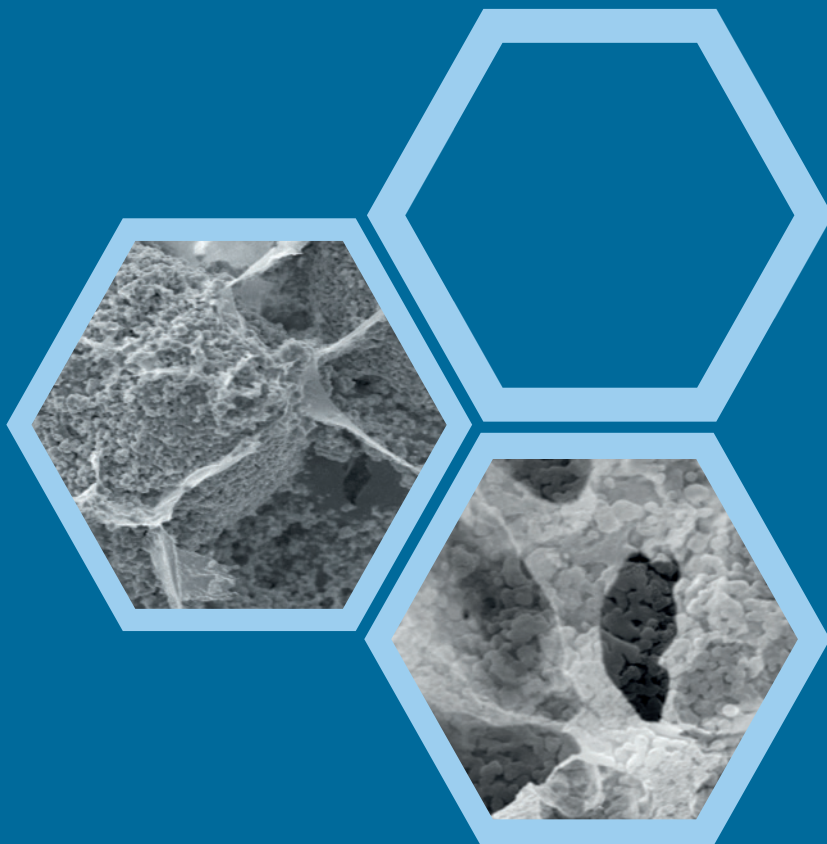٥٨ \_\_\_\_\_ **EZ Cure™**

٦٠ \_\_\_\_\_ **MBCP+™**

**Biomatlante**  
Biologics Solutions  
France

# IN'OSS (MBCP Putty)

## IN'OSS - MBCP Putty Bone Graft Substitute

|                                       |                                                                                                                                                      |
|---------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------|
| Synthetic with hydrogel               | Developed on MBCP Technology: Biphasic Calcium Phosphate (HA/ $\beta$ -TCP) + hydrogel                                                               |
| Developed to facilitate handling      | Keeps the original graft shape and bone volume                                                                                                       |
| Global Porosity                       | Interconnected network that enables the colonization of bone cells and biological fluid uniformly within the matrix.                                 |
| Macroporosity<br>>100 $\mu$ m average | Macropores are a network of interconnected spaces that promote the biological infiltration and cellular colonization by osteoblasts and osteoclasts. |
| Microporosity<br><10 $\mu$ m          | Micropores are the intercrystalline spaces where dissolution and recrystallisation occurs.                                                           |
| Osteoconductive                       | Provides a matrix for new bone growth                                                                                                                |
| Bioactive                             | For ionic exchange: TCP dissolution and bone crystal precipitation create newly bioactive interface with bone cells.                                 |
| Sterilisation                         | Steam                                                                                                                                                |
| Shelf life                            | 3 years                                                                                                                                              |
| Available in different volumes        | 0,5 mL - 1 mL<br>2,5 mL - 5 mL - 10 mL                                                                                                               |

## Moldable and Injectable Bone Graft Substitute

Building on Biomatlante's core MBCP Technology, In'Oss is a moldable synthetic bone graft, composed of Hydroxyapatite (HA), Beta Tricalcium Phosphate ( $\beta$ -TCP) and a hydrogel.

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## Resorbable Collagen Membrane

EZ Cure Membrane, flexible and adaptable

|                                   |                                                                                                                                                                                                                                                            |       |       |
|-----------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|-------|
| Origin                            | EZ Cure is a highly purified porcine collagen (fibers type I and III) resorbable membrane.                                                                                                                                                                 |       |       |
| Resorbable                        | The resorbable collagen membrane allows biological fluids diffusion for a controlled resorption within 24 weeks.                                                                                                                                           |       |       |
| Barrier function                  | <p>The cross-linked collagen membrane offers a controlled resorption time that ensures a barrier function for 12 weeks.</p> <p>This barrier effect prevents epithelial and connective tissue from colonizing the surgical site during healing process.</p> |       |       |
| Reversible                        | EZ Cure can be placed on either side thanks to its 3 layers structure. It is very easy to use.                                                                                                                                                             |       |       |
| Flexible and adaptable            | <p>EZ Cure has excellent handling properties.</p> <p>Flexible and drapable, it conforms perfectly to various bone geometries.</p> <p>Adaptable, it is easy to cut and shape. It can be sutured and pinned.</p>                                             |       |       |
| Adhesive                          | <p>Thanks to its surface morphology, the membrane presents a great tissue adherence.</p> <p>It is easy to place.</p>                                                                                                                                       |       |       |
| Sterilisation                     | Irradiation                                                                                                                                                                                                                                                |       |       |
| Shelf life                        | 3 years                                                                                                                                                                                                                                                    |       |       |
| Available in different sizes (mm) | 15x25                                                                                                                                                                                                                                                      | 20x30 | 30x40 |

### Dental Surgery:

EZ Cure Membrane is intended for use in Guided Tissue Regeneration (GTR) or Guided Bone Regeneration (GBR) procedures.

Its barrier function prevents colonization of the bone defect by epithelial and connective tissue.

It is indicated for treatment of oral and maxillofacial bone defects and periodontal diseases.

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
Optimal wound healing

Cell attachment promotion

Crosslinked for long-lasting barrier

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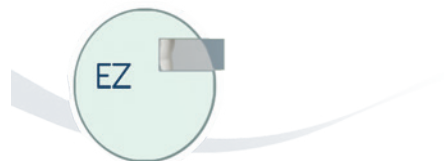## MBCP Synthetic Bone Graft Substitute



|                                     |                                                                                                                                                   |
|-------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------|
| Synthetic                           | Biphasic Calcium Phosphate (HA/β-TCP)                                                                                                             |
| Two ratio available                 | 60% HA/ 40% β-TCP<br>20% HA/ 80% β-TCP*                                                                                                           |
| Global Porosity of 70%              | Interconnected network of macropores and micropores that enables the colonization of bone cells and biological fluid uniformly within the matrix. |
| Macroporosity<br>300-600 μm average | Interconnected spaces that promote the biological infiltration and cellular colonization by osteoblasts and osteoclasts.                          |
| Microporosity<br><10 μm             | Micropores are the intercrystalline spaces where dissolution and recrystallisation occurs.                                                        |
| Osteoconductive                     | Provides a matrix for new bone growth                                                                                                             |
| Bioactive                           | For ionic exchange: TCP dissolution and bone crystal precipitation create newly bioactive interface with bone cells.                              |
| Sterilisation                       | Irradiation                                                                                                                                       |
| Shelf life                          | 5 years                                                                                                                                           |
| Available in different shapes       | granules, sticks, wedges, cylinders, blocks, custom-made shapes on request                                                                        |

## Synthetic Resorbable Bone Graft

# MBCP Technology Range

## Synthetic Bone Graft Substitute

MBCP is a bioactive calcium phosphate ceramic composed of hydroxyapatite (HA) and Beta-Tricalcium Phosphate (β-TCP).
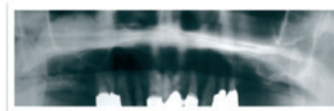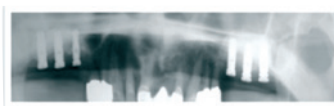
## COLLAGENE AT® Spain

Resorbable membrane for  
GUIDED TISSUE REGENERATION

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
## Properties

- Barrier effect being resorbed within **180** days
- Adaptability by adhesion to the surrounding tissue
- **Surface micro roughness** that allow the stabilization of the membrane **without suture**
- The optimal **barrier effect** is achieved in association with any grafting materials
- The membrane is sterilized by means of gamma radiation and packaged
- In a double envelope, guaranteeing a safe and protected sterilization

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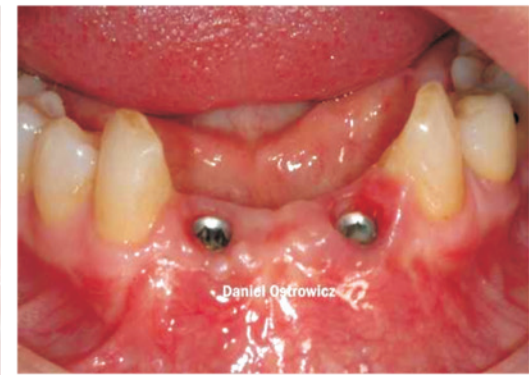### Camera Auto Focus:

Lens : Zoom 60X  
Considerable Depth of Field  
System : NTSC / PAL  
Number of Pixels : 1080 (1080 \* 1920)  
Horizontal Field of the Vision : 4.8 to 2.8  
Low Need of Illumination : 0.7 Lux  
White Balance : Automatic , Indoor , Outdoor  
Video Output : HDMI - YPBPR  
Power Waste : 6-12v DC , 3.3v

### Light:

White Natural Day Light : 6300 k - 6500 k  
Good Reproduction of Colors : 80%  
Adjustable Light Intensity until 55/000 Luxes  
No Shadow Carried  
Completely Light and Water Proof for on Optimal Asepsis  
Absence of UV Radiation and no Extra heat Contribution  
Whitout Fan and Noiseless  
Source : 2LED \* 3.15 w , Consumption 9.45 w  
Life Expetancy of 50/000 Houres  
Energy Supply 110-220v , Autonomous  
Stabilized Converter

### Whole Set:

Dimension of the Head : 272 \* 80 \* 146mm  
Scale of the Arm : 1470mm  
Weight of the Head : 1.75 kg  
Weight of the Arm : 6.90 kg  
Diameter of Axis : 35mm  
Color : RAL9016  
The set is Delivered With Arms  
Warranty 1Year

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## QuickWhite In-Surgery/Laser (Office) ( 3 patients kit )

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Our In-surgery kit comes as a kit form and conforms to the UK & European teeth bleaching regulations.

The famous and best selling in-surgery/Laser whitening kit for teeth Bleaching.

Contains 3 pots or syringes of QuickWhite for 3 Patients, 2 Paint-on dam 3g syringes, brushes, pumice, vitamin Es and applicators and colour instruction.

Laser kit QW1001L have Potassium Nitrate 10% wv. Pots for use with 6% - 38% hydrogen peroxide.

kits with higher peroxide UK/EU allowed limits are for overseas. Dental use only.

## Carbamide Peroxide Syringes (home) ( 1 patient )

Comes in either 5 syringes or 10 syringes packs

**Night Time:** 10% Carbamide kit, large 3cc syringes packs with re-sealable tips, patient's whitening guide and whitening instructions.

**Night/Day time:** 16% Carbamide kit, large 3cc syringes packs with re-sealable tips, patient's whitening guide and whitening instructions.

22% Carbamide kit, large 3cc syringes packs with re-sealable tips.

30% Carbamide kit, large 3cc syringes packs with re-sealable tips.

kits with higher peroxide UK/EU allowed limits are for overseas.

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
**MULTIFUNCTIONAL PEDAL CONTROL (optional)**


with reverse direction. on-off pump, start speed regulation and possibility to change functions P1/P2/P3/P4. The pedal controls Micromotor start. Suitable for operating theatres and sterile environments.

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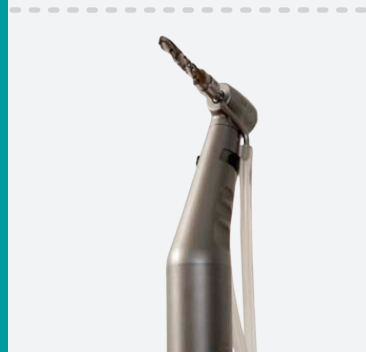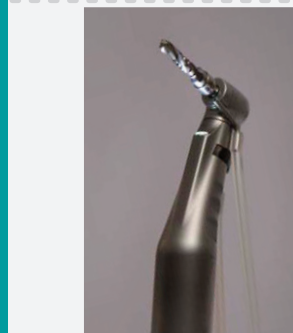## AREA OF USE:

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- ✓ Sinus Lift with crestal approach
- ✓ Extraction
- ✓ Apicectomy
- ✓ Cystectomy
- ✓ Bone Harvesting
- ✓ Ridge Expansion
- ✓ Sinus Lift
- ✓ Osteoplasmy, Osteotomy
- ✓ Implant site preparation
- ✓ Tooth preparation
- ✓ Orthodontic Surgery

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High power handpiece  
 with internal spray.

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Pump with easy  
 insertion of the  
 irrigation set.

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Kit with 6 insert tips  
 (5 sharp, 1 smoothing).

# SURGYBONE TIPS

A new generation in ultrasound surgery

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
RECOMMENDED SETTING:

PWR 30  
VIBRA 100

• If you exceed the recommended settings, we disclaim responsibility and guarantee

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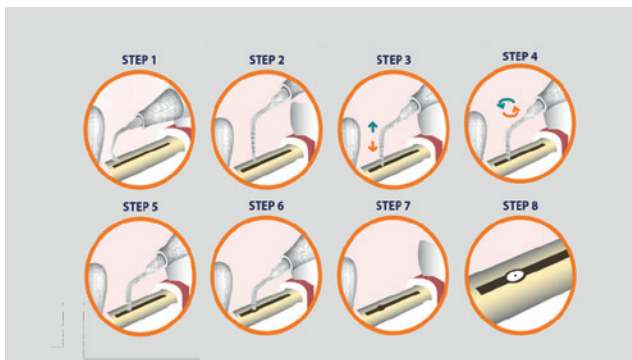## CGF KIT

#### AREA OF USE


- ✓ Orthopedic surgery;
- ✓ Maxillo facial surgery;
- ✓ Oral surgery;
- ✓ Ophthalmology;
- ✓ Cosmetic surgery;
- ✓ Sport medicine;
- ✓ Dermatology;
- ✓ Gynecology;
- ✓ Neurosurgery.

## ROUND UP

Patent pending - Integrated Mixing

- Medical device for intrinsic and extrinsic molecular blend, mixing with altering the geometric dimensions of autologous, heterologous or synthetic materials, for medical use;
- The system is used with liquid, semi-liquid or solid materials;
- Perfect, homogenous mixing free from atmospheric contamination;
- Automatic empty in just a few second: max 16 seconds.
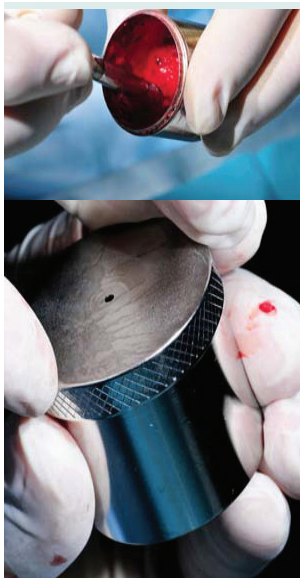
### CGF - Concentrated growth factors -

It represents a new generation of platelet able to hold inside a higher concentration of autologous growth factor. It is isolated from blood samples through a simple standardized separation protocol, which is performed by means of a specific centrifuge without the addition of exogenous substances.

The patent kit includes the necessary for the blood collection and specific dappens to facilitate the processing of hematic components without any contamination risks.

Then, the kit is equipped with specific instruments to produce a fibrin membrane and to insert the clot inside the implant site.

### • POWER SOURCE

✓ 100: 115+/-10% | 50/60 Hz

### • WEIGHT

✓ 19 Kg

### • DIMENSIONS

✓ 320 P X 320 L X 300 H

### • NOMINAL POWER CONSUMPTION

✓ 320 VA

### • RPM

✓ 3600 r.p.m

### • TIMER

✓ 1-16 SEC.

### THECNICAL DATA

- Equipped with a decontamination cycle with UVC reflected light;
- Cycle duration: 5 minutes.

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